



CE



Mewa

Operation and maintenance manual

Wydanie 2 z dn 2020-12-18

A modern stroller of excellent driveability MEWA



ATTENTION! Users with sight, reading or cognitive impairment are requested to report to the Manufacturer in order to obtain information for the user in an appropriate format.

Purpose of the Manual

The operation and maintenance manual contains basic information necessary for the correct use of the Mewa stroller, including the principles of safe use, description of its construction, operation and adjustment, preparation of the stroller for operation, how to maintain and clean the stroller, and warranty conditions.

It is designed for persons who directly care of children with disabilities (assistants child - parents, guardians), which will be transported in a stroller. It is also intended for doctors and physiotherapists operating the stroller.

ATTENTION! Keep this manual for future reference!

Keep the manual in an easily accessible place.

Manufacturer's contact details:

- Postal address: MyWam sp.j., ul. Szczecińska 10, 41-516 Chorzów, Poland
- contact phone number: +48 32 733 11 31
- e-mail address: mywam@mywam.pl



ATTENTION! Before using the stroller, each User is obliged to read and follow this manual.
ATTENTION! This manual is an integral part of the product and must always be attached to the product being sold.

User and maintenance manual version: EN

Edition: II 2020-12-18

MyWamKupiec, Bartold, Angres sp.j. reserves the right to make technical and commercial changes to the content of the Operation and Maintenance Manual without prior notice. Each change will be determined by the date the instructions were last updated .

Table of contents

▪ Introduction.....	5
1. Intended use of the Mewa stroller	5
2. Basic information about the Mewa stroller.....	6
3. Quick start.....	6
▪ Safety.....	7
4. Important instructions and safety rules.....	7
5. Nameplate.....	13
6. Warning labels	13
7. Symbols designation	13
▪ Product	14
8. General construction of the stroller	14
9. Configuration and use – stroller frame	15
9.1. Frame folding mechanism	15
9.2. Unfolding the stroller frame. Transport lock.....	15
9.3. Folding the frame.....	15
9.4. Swivel front wheels, solid. Locking the direction of travel of the front wheels.....	16
9.5. Rear wheels, solid. Installation and removal of rear wheels on the axle.....	16
9.6. Rear wheel amortization system	17
9.7. Ergonomic handle for pushing the stroller. Mechanism and angle adjustment of the handle	17
9.8. Central brake. Applying and deactivating the brake	18
9.9. Tilting foot.....	18
9.10. The transport handle	18
9.11. Footrest. Adjustment and use.	18
10. Configuration and use - therapeutic seat with backrest.....	19
10.1. Assembly and disassembly of the seat to the frame	19
10.2. Adjustable seat, „reclining back”	19
10.4. Adjustable backrest.....	20
11. Basic equipment	20
11.1. Five-point belts	20
11.2. Safety barrier.....	21
11.3. Adjustable canopy with window	21
11.4. Shopping basket.....	21
11.5. Abduction and stabilizing belts	22
11.6. Wedge	22
11.7. Lininig pads.....	22
11.8. Upholstery.....	22
12. Additional equipment	23
12.1. Adjustable head supports - ↔ and ⇅	23
12.2. Adjustable chest and loin pads - ↔ and ⇅	23
12.3. Therapeutic table	23
12.4. Seat minimizing kit.....	24
12.5. Travel bag.....	24
12.6. Sun umbrella	24

12.7.	Winter gloves	24
12.8.	Winter sleeping bag	25
12.9.	Mosquito net.....	25
12.10.	Raincover.....	25
13.	Technical specification of the Mewa stroller	26
13.1.	User parameters	26
13.2.	Material composition of the stroller.....	26
14.	Delivery and tools.....	26
17.	Using the stroller	28
17.1.	Carrying the stroller	28
17.2.	Getting on and off the stroller	28
17.3.	Correct position in a stroller	29
17.4.	Riding in a stroller	29
17.5.	Moving on slopes	29
17.6.	Overcoming thresholds or curbs.....	29
17.7.	Transport in a car and other vehicles	29
17.8.	Cleaning and care.....	29
18.	Storage of the stroller	30
18.1.	Preparation of the stroller for long-term storage.....	30
18.2.	Corrosion protection.....	30
19.	Re-use.....	31
23.	Notes	34
24.	List of authorized service points	35
25.	Warranty Card	36
26.	Warranty conditions.....	36

▪ Introduction

Thank you for purchasing special **Mewa** stroller for young people. We hope that it will be the cause of many wonderful moments with your child and will be helpful every day. We make every effort to ensure that our products meet your expectations and the highest standards. We encourage you to contact our employees and follow us on our website www.mywam.eu and on social media. There you can find out many interesting things about our products.

1. Intended use of the Mewa stroller

- The Mewa Special Stroller is designed to transport a child with neurological, orthopedic or other conditions that hinder or prevent him from moving independently.
- The specialistic Stroller Mewa was designed for children who suffer from neurological and orthopedic problems, which make it hard or impossible for them to move independently. The stroller is equipped with adjustments which make it possible to safely and comfortably adapt it to the child's individual needs.
- The Mewa stroller has been designed to carry one child.
- The Mewa stroller, through the use of many regulations, allows comfortable and safe adaptation to the individual needs of the child.
- **The Mewa stroller is designed to carry a child whose weight does not exceed 35kg, height from 90cm to 125cm.** Statistically, they are children from 3 to 8 years old.
- When choosing a stroller, however, the most important are the individual physical and mental characteristics of the child, its dimensions and weight, housing conditions, environment and, above all, the individual's illness.
- Among the disease entities of a particular child, among others, there can be mentioned such diseases as: children cerebral palsy, meningitis, condition after spinal injury, condition after craniocerebral trauma, condition after stroke, paralysis of various origin and other dysfunctions damaging the locomotor and preventing independent movement of the child.
- Use the stroller as directed by a specialist.
- The Mewa stroller is intended for both indoor (e.g. apartment, rehabilitation center) and outdoor (e.g. park, pavement), on hard and level ground, on which all four wheels touch the ground and their contact with the ground is sufficient to properly drive the wheels.
- The Mewa stroller is intended to be controlled and propelled only by an assistant (parent, guardian).
- The user of the Mewa stroller is both the child who will be transported in it and the child's assistant (his parent, guardian, doctor, physiotherapist) who will transport the child in the stroller and look after him.
- Mewa stroller with proper care and maintenance maintains its quality for a period of 3 years and this is the expected period of its use.
- The risks associated with using the stroller are described in detail in the chapter: "Important instructions and safety rules".


2. Basic information about the Mewa stroller

The Mewa youth stroller has many options, e.g.:

- backrest angle adjustment (see 10.3)
- backrest and seat angle adjustment "Cradle" (see 10.2)
- seat depth adjustment- seat minimizing kit (see 12.4)
- footrest adjustment (see 9.11)
- torso stabilization (see 12.1 & 12.2)
- optional equipment (see 12)

3. Quick start



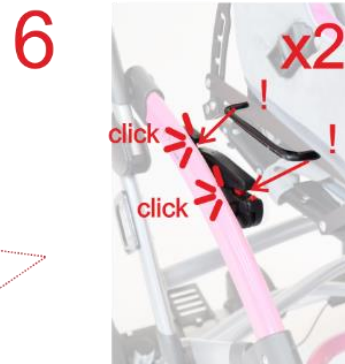






To quickly use the stroller it must be unpacked and prepared for use according to item 15 Assembly and disassembly of the stroller and item 16 Preparation for use. The following is a pictorial diagram, but it will not replace careful reading of the complete manual.



▪ Safety




4. Important instructions and safety rules

General tips:


-  • **WARNING!** Risk of accident - before using the stroller, the user is obliged to read and follow the entire contents of this manual.
-  • **WARNING!** Risk of accident - Before starting and while using the stroller, it is essential to comply with the safety instructions and rules described in this manual, and comply with the other instructions described in this manual.
-  • **WARNING!** Failure to follow the manual, in particular those regarding the moving parts of the stroller, may be dangerous for the user.
 - Always read the operating and maintenance manual before using the stroller.
-  • **WARNING!** Never leave your child unattended in a stroller, as there is a risk of injury and there is a risk of accident!
-  • **WARNING!** The product contains small parts that can be dangerous if swallowed by a child!
-  • **WARNING!** Risk of injury - always use all available security systems for the child using a stroller. These include: five-point belts, seat minimizing kit, abduction and stabilizing belts, wedge, safety barrier or in its place a therapeutic table, footrest and straps for placing feet on the footrest.
-  • **WARNING!** There is a risk of burns when using the stroller - be careful when transporting a child in a stroller at high and low temperatures and in contact with the skin, e.g. in the sun, in cold weather, in a sauna, because the temperature of the stroller surface may rise when the stroller is exposed to solar radiation or the stroller surface temperature may drop when exposed stroller for frost.
-  • **WARNING!** When using the stroller or operating the stroller or folding the stroller, there may be a risk to the health of users and people in the immediate vicinity, caused by the dimensions of fixed and movable gaps and openings. Be particularly careful when doing this.
-  • **WARNING!** When using and operating the stroller, as well as when folding and adjusting the mechanisms, there may be a risk of entrapment and squeezing of parts of the body of the user or accompanying person in the openings and gaps between elements.
 - Be particularly careful when doing these operations.
 - After completing these operations, stabilize the positions of the individual components of the stroller by tightening the nuts, knobs and screws thoroughly. Each time you perform the above steps, you must first pull the child out of the stroller!
 - **ATTENTION!** Never exceed the maximum load.
 - **ATTENTION!** For Mewa stroller with load the maximum total load amounts to 35 kg. This value should not be exceeded.
 - **ATTENTION!** The maximum load of the travel bag is 2 kg.
 - **ATTENTION!** The maximum load of the shopping basket is 3 kg.
 - **ATTENTION!** The maximum therapeutic table load is 3 kg.
 - **ATTENTION!** The back pocket in the adjustable booth with the window is only used to place raincover or a mosquito net.

- **ATTENTION!** Depending on the adjustments of the Mewa stroller - it concerns the position of the push handle, backrest height adjustment, backrest angle adjustment, seat depth adjustment, seat angle adjustment, footrest angle adjustment and footrest settings - this stroller may exceed the recommended maximum dimensions.
- **ATTENTION!** In the event of any medical incident related to the presence of a Mewa stroller, you must immediately notify the Manufacturer in writing by contacting via registered mail or email with confirmation of receipt using the addresses below:
Postal address: MyWam sp.j., ul. Szczecińska 10, 41-516 Chorzów, Poland
E-mail address: mywam@mywam.pl
- The therapy seat is not suitable for children under 3 years old
- It is forbidden to use the stroller as a ladder or as a base for climbing onto it to reach higher.
- The stroller must not be used in a way for which it has not been designed.

Tips in particular on construction and regulation:

-  • **WARNING!** All adjustments of the stroller and its elements, all settings, changes, fastenings and other activities related to matching the stroller should be carried out **WITHOUT A CHILD** placed in the stroller. Failure to do so may cause serious injury to a child and other people!
-  • **WARNING!** Adjustments made in the stroller, made incorrectly or not in accordance with the instructions in this manual, pose a serious danger to the child and other persons!
-  • **WARNING!** Due to the large adjustment possibilities of the stroller, when adjusting the seat angle and angle adjustment of the backrest, it is possible to set the position in which the child's head will be below the pelvis line. It is forbidden to use this position, as this could lead to the risk of aspiration or choking of the child.
- **ATTENTION!** Any change in the seat tilt changes the overall centre of gravity, which affects the stability of the stroller.
- **ATTENTION!** Every change in the backrest's tilt causes a change in the overall centre of gravity, which affects the stability of the stroller. Any change in backrest tilt requires five-point belt re-adjustment!
- **ATTENTION!** The central brake is not used to slow down the stroller during movement and should only be used to immobilize the stroller.
- The front steering wheels, if the direction of travel lock is off, may vibrate and lead to a sudden stop of the stroller and its falling forward when the surface is uneven. Therefore, with an uneven surface, it is recommended to block the front direction wheel lock and push the truck at a reduced speed.
- All adjustments and settings to be performed under the supervision and care of a doctor and a physiotherapist.

Information particularly regarding assembly and disassembly:

-  • **WARNING!** To avoid injury, make sure the child is at a far and safe distance when the stroller is assembled and disassembled.

- ⚠ • **WARNING!** The stroller must be fully assembled before use. In case of any problems at this stage, please contact the Seller and the Manufacturer urgently.
- ⚠ • **WARNING!** When mounting the stroller and when mounting the seat to the stroller frame, pay special attention to mount the seat in the designated place of the frame and secure all belts.
 - Fastening the seat to other elements of the stroller frame is incorrect and may cause an accident!
- ⚠ • **WARNING!** It is forbidden to place additional weight on the stroller, which in particular applies to placing additional weight on the back of the back e.g. backpack or bag. This may cause the Mewa stroller to tip over and cause user injury! Does not apply to original additional accessories dedicated to the Mewa stroller not exceeding the permissible load.
- ⚠ • **WARNING!** When assembling and disassembling the stroller, keep a safe distance from moving parts of the stroller, as there is a risk of fingers, hands or other body parts being trapped.
- ⚠ • **WARNING!** When assembling and disassembling the wheels to the stroller frame, make sure and check that they are correctly fitted, as there is a risk of injury.
 - **ATTENTION!** Do not unfold the stroller frame if the transport lock is in place, as it may damage it.
 - All additional elements attached to the stroller affect its stability.

Tips in particular regarding preparation for use:

- ⚠ • **WARNING!** Before using the stroller, make sure that the stroller is properly assembled and all mechanisms, catches and fasteners are properly assembled, closed and prepared for use. In particular, this applies to: the mechanism for folding the frame, the mechanism for attaching the seat to the frame of the stroller and the installation of all four wheels.
- ⚠ • **WARNING!** Before riding the Mewa stroller, make sure that the child is secured and correctly adjusted and fitted with five-point belts, as there is a risk of injury.
- ⚠ • **WARNING!** The frame folding mechanism is a very important element of the construction of the stroller, so each time before using the stroller, pay special attention to the absence of cracks and damage to the folding mechanism and its proper operation. In the event of any suspicion you must not use the stroller, as this may endanger the user, child's assistant and other people! In this situation, it is recommended to urgently contact the Seller and the Manufacturer for advice.
- ⚠ • **WARNING!** Correct use of the frame folding mechanism is important for the safety of the child and other people.
 - **ATTENTION!** The upholstery must always be dry before using the stroller.
 - **ATTENTION!** Before each use of the stroller, the activities listed below must be carried out, and in the event of irregularities the stroller must not be used. In this situation, it is recommended to urgently contact the Seller and the Manufacturer for advice.



The operations before each use of the stroller include:

- Check the stroller for visible damage or other defects that could endanger the child during stroller operation.
 - Check that the stroller frame has been unfolded correctly and the frame folding mechanism has been correctly closed (locked).
 - Check that the seat has been correctly mounted to the stroller frame and that the seat attachment mechanism to the stroller frame has been properly closed (locked).
 - Check the correct assembly of all four wheels.
 - Check the good condition of the tread tires on all four wheels.
 - Check good central brake operation.
 - Correctly make all adjustments and adjust the stroller to the child's needs, including setting the appropriate depth of the therapeutic seat and height of the backrest, appropriate inclination of the therapeutic seat and inclination of the backrest, set the appropriate inclination of the footrest, adjust the appropriate length of five-point belts.
 - Check that all adjustments are properly seated and secured.
 - Check that the stroller upholstery is dry.
 - Secure the child with five-point belts to prevent the child from falling out and sliding out of the stroller.
 - Secure the child with all available stabilizing elements, such as five-point belts, seat minimizing kit, pelvic support belts, abduction and stabilization belts, wedge, footrest, footrest straps, safety barrier or therapeutic table. This is to protect the child from falling out and sliding out of the stroller and to protect the child's feet from falling from the footrest.
- **ATTENTION!** The stroller must not be used if the technical condition of the stroller or a single element of the stroller is bad or if the user has any doubts about it, the stroller is damaged or there is only a suspicion that damage has occurred or the stroller is not working properly. In this situation, it is recommended to urgently contact the Seller and the Manufacturer for advice.




Information particularly regarding the handling of the stroller:

- **ATTENTION!** It is forbidden to carry or lift the stroller with the child inside.
- **ATTENTION!** Due to the weight of the stroller and the inability to detach the seat from the frame, lifting the stroller is discouraged carrying it "in hands". If necessary, remove all accessories mounted additionally to the stroller and the rear wheels. The stroller should be lifted and carried by 2 people. Alternatively, you can lock the front wheels and pull the stroller behind you instead of carrying it.
- **ATTENTION!** Do not lift the stroller by the footrest or safety barrier, therapy table or any other additional accessories.




Information in particular regarding getting in and out of the stroller and the correct position in the stroller:

-  • **WARNING!** Do not stand on the footrest, as there is a risk of the stroller tipping over, injury to the user and risk of damage to the product.
-  • **WARNING!** Make sure that when getting in or out of the stroller, the child and no other person stand on the footrest, as this may cause damage and injury to the user!
- **ATTENTION!** When placing and removing a child from the stroller, the central brake must be activated.


Tips in particular regarding stroller driving:

-  • **WARNING!** The stroller must not be released by itself or pushed away, as it may cause an accident.
-  • **WARNING!** This product is not suitable for running or skating.
-  • **WARNING!** Use of the product for sport purposes is prohibited.
- **ATTENTION!** Make sure that during use of the stroller, the child's limbs do not get between the wheel spokes and the adjustable elements of the stroller.
- **ATTENTION!** When passing through a door or similar narrowing limiting the direction of the stroller, make sure there is adequate space on both sides of the stroller, which prevents your hand or other body parts from being trapped, as well as preventing damage to the stroller.
- When traveling on public roads, follow the Highway Code.

Tips in particular on moving on slopes:

-  • **WARNING!** It is forbidden to leave a child in the stroller without supervision of an adult and responsible assistant (parent, guardian), which in particular concerns leaving the child in the stroller on a slope, as it may cause an accident and injury to the user!
-  • **WARNING!** When using ramps, do not exceed the maximum load of the ramps, as this may cause injury to users and people around you.
-  • **WARNING!** When using ramps, make sure that the ramp is suitable, that the wheels of the stroller are high enough to overcome the driveway, that the stroller frame does not touch the driveway, and that help should be provided to another person as there is a risk of the stroller falling over, damage to the product and injury to the user.
- **ATTENTION!** Always remember that the location of the stroller's centre of gravity changes when overcoming unevenness or slopes. In this case, be extremely cautious and secure the stroller.

Tips in particular on overcoming thresholds or curbs:

-  • **WARNING!** It is forbidden to climb or slide the stroller down the stairs, both with the child and the stroller itself, because it may cause an accident and injury to the user!

Information particularly regarding the transport of the stroller in a car and other vehicles:

- **WARNING!** It is forbidden to transport a child in a stroller and at the same time in any other means of transport, e.g. in a car, on a train, in a tram, on a bus, on a plane, on a ship or other, because it may cause injury.
- **WARNING!** It is forbidden to use the stroller seat as a car seat, as this may cause injury.
- **WARNING!** Risk of accident and injury to the user and other people - when transporting the stroller in other means of transport, e.g. in a car, bus, train, make sure that the stroller is correctly fastened to prevent injury to passengers during collision or sudden braking.
 - **ATTENTION!** When transporting the stroller in other means of transport, e.g. in a car, on a bus, on a train, the stroller must be transported without a child and the stroller must be well protected against any damage.
 - **ATTENTION!** The manufacturer is not liable for damages resulting from inadequate product protection during transport.
 - **ATTENTION!** When transporting the stroller in other means of transport, e.g. in a car, on a bus, on a train, do not use one seat belt to secure both the passenger and the stroller at the same time.
 - Use the transport lock when transporting the stroller and storing it.

Tips especially on maintenance and service:

- **ATTENTION!** The Manufacturer is liable for defects or technical damage to the stroller only caused by proper use of the stroller. The responsibility for the proper functioning of the stroller rests with the Manufacturer only if you strictly follow all the Manufacturer's instructions.
- **ATTENTION!** The Manufacturer is responsible for the proper functioning of the stroller only if original parts and original additional accessories are used.
- Only original spare parts that are supplied or recommended by the Manufacturer may be used in the stroller.
- Repairs to the stroller carried out by persons other than the authorized Service (see Chapter 19.) or the Manufacturer will void the warranty. The Manufacturer is not responsible for personal injury or damage to the product that may result from repairs carried out by persons other than the authorized Service centre (see Chapter 19.) or the Manufacturer.
- The stroller should be cleaned and looked after regularly. Do not use abrasive cleaning agents.
- The Manufacturer is responsible for the proper functioning of the stroller only if it is purchased in a specialized medical store.
- The Manufacturer is not responsible for the appropriate selection of accessories for the given user.

5. Nameplate

- The nameplate contains basic information about the product: name, model, serial number (identification), maximum load, date of manufacture, confirmation of product compliance for medical devices, warning about proceeding in accordance with the Operation and Maintenance Manual and manufacturer's data (see Fig. 5.1).
- The nameplate is located on the crossbeam at the rear of the stroller frame, on its left side, glued from above.

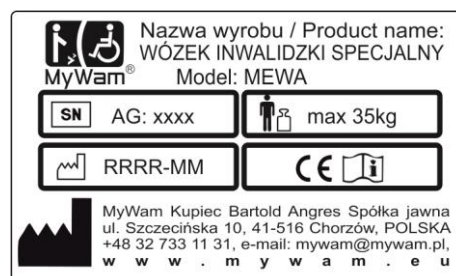


Fig 5.1

6. Warning labels

- The warning label contains important warnings for the user (see Fig. 6.1).
- The warning label is located on the rear/side part of the frame on both sides.

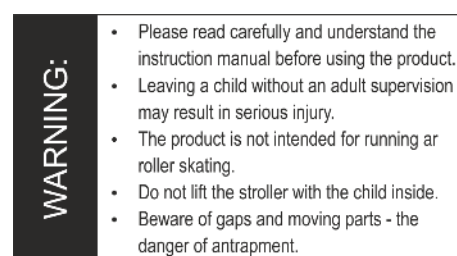


Fig 6.1

7. Symbols designation

	Product serial number		Nazwa producenta
	Permitted total user weight		DESCRIPTION " WARNING! " and the symbol in this manual are used to enhance the reader's attention to the content marked with this sign. Non-compliance with content marked with this description may endanger the health and safety of the user and accompanying persons, and may also damage the product.
	Date of manufacture		
	Conformity mark according to the Medical Devices Directive No. 93/42 EEC Annex VII		
	Warning about proceeding in accordance with the product Operation and Maintenance Manual		

▪ **Product**

8. General construction of the stroller

The Mewa Special Stroller for children consists of two main parts: the frame of the stroller with seat and backrest and very rich additional equipment that should be selected according to the individual needs of users.

• **Standard equipment**

The basic elements of the stroller frame include:



- 1- Adjustable canopy with window - see 11.3
 - 2- Ergonomic handle for pushing the stroller - see 9.7
 - 3- Five-point belts - see 11.1
 - 4- Backrest and seat angle adjustment mechanism - see 10.2
 - 5- Transport lock - see 9.2
 - 6- Adjustable footrest with straps - see 9.11
 - 7- Rear wheel with amortization system - see 9.5
 - 8- Lock of the direction of travel of the front wheel - see 9.4
 - 9- Front steerable wheel, solid - see 9.4
 - 10- Safety barrier - see 11.2
 - 11- Wedge - see 11.6
 - 12- Backrest and seat upholstery - see 17.8
 - 13- Central brake - see 9.8
 - 14- Tilting foot - see 9.9
 - 15- Rear wheel amortization system - see 9.6
 - 16- Lininig pads - see 11.7
- *) option

• **Optional equipment**

Its elements may include:

- Therapy table - see 12.3
- Abduction and stabilizing belts- see 11.5
- seat minimizing kit - see 12.4
- Adjustable chest and loin pads - see 12.1&12.2
- Travel bag - see 12.5
- Shopping basket - see 11.4
- Raincover - see 12.10
- Winter gloves - see 12.7
- Winter sleeping bag - see 12.8
- Mosquito net - see 12.9
- Sun umbrella - see 12.6



9. Configuration and use – stroller frame

9.1. Frame folding mechanism

The frame folding mechanism is located inside the frame and on the frame in two places in its right part and in three places in its left part. Describing these places from above, from the side of the handle to push the stroller, these are:

- right and left bend - an element that allows folding the upper frame to reduce its dimensions
- mechanism of folding the frame with a lock



Fig 9.1

Places of folding/unfolding

9.2. Unfolding the stroller frame. Transport lock.

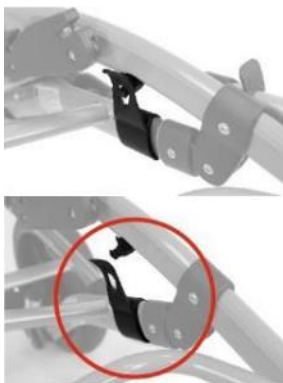


Fig. 9.2.1

Transport lock



Fig. 9.2.2

Unfolding the

To unfold (open) the stroller frame: lay the stroller frame horizontally on the ground and pull off the transport lock (Fig. 9.2.1). Then gently pull the backrest to the upright position (Fig. 9.2.2.). To set the stroller, pull the handle carefully by lifting it slightly up while bringing the stroller to the desired position until the frame folding lock engages. The frame of the stroller will not unfold properly until the frame folding mechanism has been properly closed. This is evidenced by the clear "click" sound coming from the well-closed mechanism of folding the frame

9.3. Folding the frame

To fold (close) the stroller frame: push the small, red button (Fig 9.3.1) and release the frame folding mechanism lock by pulling it up and pull both frame folding handles up (Fig 9.3.2) until the frame is released. Then fold the stroller frame folding it forward in the direction of the front wheels. Make sure that the push handle and footrest are in a non-colliding position. Then pull the backrest release mechanism. Put on a transport lock.



Fig. 9.3.1

Red button

Fig. 9.3.2

Folding frame mechanism

9.4. Swivel front wheels, solid. Locking the direction of travel of the front wheels

The used, durable wheels guarantee the highest quality of stroller driving, driving pleasure and elegant appearance. The tires are tubeless, filled with a special gel that has shock absorbing properties very similar to those of inflatable wheels. The wheels are factory-fitted to the front wheel forks, which are finished with top steering axles from above.

Both front wheels can be set for straight-ahead or unlocked for turning, which makes it much easier to maneuver the stroller. Thanks to the locking of the direction of travel of the front wheels and the possibility of putting it on or taking off, there is clearly greater comfort for both the child and the assistant when driving the stroller.

It is recommended to remove the lock and allow the front wheels to turn while driving the stroller on a completely smooth surface, e.g. in shopping centers, school.

- **Assembly.** To attach the wheels, insert the wheel axle into the folding system slot until you hear a 'click'. The audible clicking sound is the indication that the assembly has been done properly. For dismantling the front wheels press the button on the front wheel folding system (Fig. 9.4.1) and eject the wheel axle from the front wheel folding system slot.



Fig 9.4.1
Quick release button

For dismantling the front wheels press the button on the front wheel folding system (Fig. 9.4.1) and eject the wheel axle from the front wheel folding system slot.

- **Direction lock.** In the event of uneven surfaces, moving on sloping surfaces, when crossing thresholds or curbs, it is recommended to install the lock and block the turning of the front wheels. To block the front wheels, turn the lock knob left. Turn the lock knob right, to unlock. (Fig.9.4.2)



Fig 9.4.2
Direction lock knob

9.5. Rear wheels, solid. Installation and removal of rear wheels on the axle

The used, durable wheels guarantee the highest quality of stroller driving, driving pleasure and elegant appearance. The tires are tubeless, filled with a special gel that has shock absorbing properties very similar to those of inflatable wheels.

Thanks to the possibility of dismounting both the wheels and their axle, the trolley is smaller after folding, which facilitates its transport and storage.

The wheel assembly mechanism is mounted inside the rear wheels and ensures quick assembly and disassembly of the rear wheels. It is started by pressing the quick assembly clamp (Fig. 9.6.1) towards the center of the wheel, after which the rear wheel can be mounted on the rear wheel axle or the rear wheel can be removed from the rear wheel axle.



Fig 9.5.1
Quick mount buckle

- **Assembly.** When installing the rear wheel on the rear wheel axle, press the clamp in the center of the wheel and put the wheel on the beginning of its axle. Then stop pressing the clamp and continue

sliding the wheel onto its axle until it clicks. The correct wheel placement is demonstrated by the clear "click" sound coming from the quick axle rear axle mounting mechanism.

- Dismantling the wheels from the axle. In order to remove the rear wheel from the rear wheel axle, press the clamp in the middle of the wheel and slide the wheel out of its axle. – Fig 9.5.2.



Fig 9.5.2
Wheel disassemble

9.6. Rear wheel amortization system

The stroller is equipped with spring-loaded cushioning of the frame and rear wheels, located at the ends of the rear legs of the frame, when in contact with the crossbeam. Amortization works automatically. The adjustment is done by activating the up-down lever located from the bottom of the shock absorber. In order to facilitate the switchover of settings by the lever, simultaneous pressure on the shock absorber spring is allowed. When the lever is shifted, the shock absorber spring will be compressed (option for hard ground) or extended (soft ground).

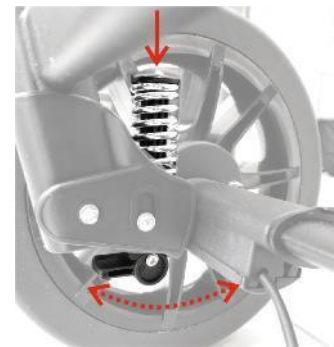


Fig 9.6
Wheel amortization system

9.7. Ergonomic handle for pushing the stroller. Mechanism and angle adjustment of the handle

The handle is used to push the stroller, providing a firm grip in various driving conditions. It has the ability to adjust the height of the handle for the child's assistant who drives the stroller - Fig. 9.7

- Adjustment mechanism. These are two buttons on the outside of the frame, located on the joints of the handle to push the stroller.
- Angle adjustment of the stroller push handle. To change the angle of the stroller push handle, push both buttons of the angle adjustment mechanism of the push handle and change its angle. After setting the handle in the correct position, stop pressing the buttons, which will automatically lock the push handle in the desired position.



Fig 9.7
Handle angular adjustment

9.8. Central brake. Applying and deactivating the brake

The central brake pedal (Fig. 9.8), which when pressed blocks the two rear wheels of the stroller directly, is located in the middle of the crossbeam, at the rear of the stroller frame.

- Central brake on and off

To apply the central brake, press the brake pedal down towards the ground until it clicks. The correct activation of the brake is indicated by a clear "click" sound from the central brake.

To disengage the central brake, pull the brake pedal up.



Fig 9.8
Central brake

9.9. Tilting foot

The tilting foot (Fig. 9.9) is located on the right side of the crossbeam, at the rear of the stroller frame. In order to use it, press the tilting foot and at the same time pull the handle to push the stroller, which in turn will cause the front of the stroller to rise up together with the front wheels. Using the tilting foot makes it easier for the stroller to climb onto small steps or curbs.



Fig 9.9
Tilting foot

9.10. The transport handle

The transport handle is located under the shopping basket. Fig. 9.10



Fig 9.10
Transport handle

9.11. Footrest. Adjustment and use.

The adjustable footrest is included in the standard equipment of the Mewa stroller.

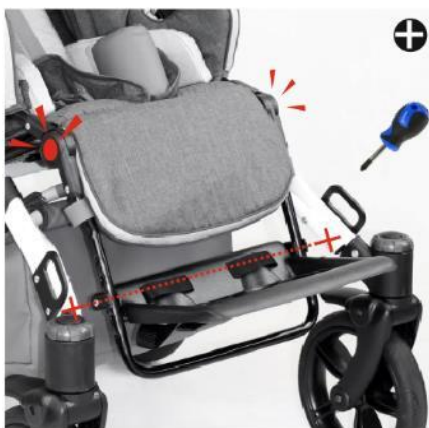


Fig 9.11
Footrest adjustment

- To change the height of the footrest, unscrew the knobs on the two screws that secure the footrest to the frame on both sides of it, then set the footrest at the appropriate height and screw it back with both screws

- To adjust the angle of the footrest, press the buttons marked in the picture, adjust the footrest and release the buttons.

- Use of a footrest. It is recommended to position the footrest at such a height that the child's feet can rest freely on it with their entire surface, while the child's thighs are parallel to the seat surface. There are straps in the footrest plate to secure the child's feet. It is recommended to put them on when using the stroller.

Normally the footrest plate is in a horizontal position. However, during transport, storage, rehabilitation, getting on and off the

Mewa stroller, the footrest plate must be folded and turned towards the footrest frame.

10. Configuration and use - therapeutic seat with backrest

10.1. Assembly and disassembly of the seat to the frame

The two frame hooks that are attached to the seat are very important parts of the stroller. They serve the purpose of attaching the seat to the stroller's frame. The proper use of the hooks is essential for the safety of the child as well as that of other people. Make sure there aren't any cracks and damage in the stroller's frame. Also, check if it works properly. Do not use the stroller in case of any suspicion as to the proper operation of the stroller. Keep the safety of your child in mind at all times.



Zdj 10.1.11

Seat mount (to frame)



Zdj 10.1.2

Seat release button

- **Seat assembly.** The seat can face the front or the back. To install the seat, place it on the stroller's frame hooks. Press until you hear a 'click' twice. Fig 10.1.1.

- **Dismantling the seat.** To dismantle the seat, press two colorful buttons situated on the lower part of the stroller's frame hooks. Next, eject the seat from the frame hooks. Fig 10.1.2.

10.2. Adjustable seat, „reclining back”

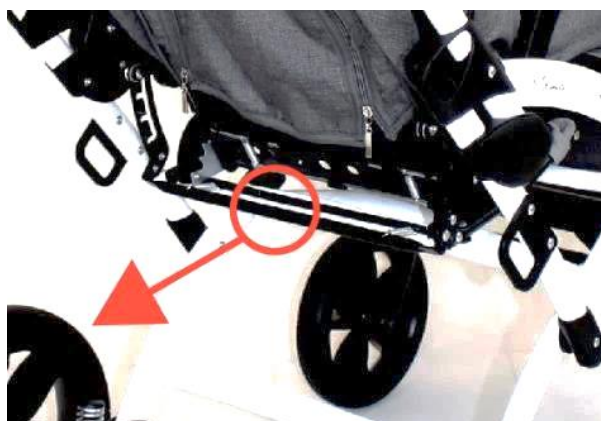


Fig 10.2

Adjustable seat

The stroller has a seat that can be adjusted to four different positions including a bassinet position. To change the position, pull back the steel handle at the back of the seat and slowly change the position of the seat. After securing the desired position, let go of the handle.

Caution! Each change of the seat position alters the general center of gravity, which impacts the stability of the stroller.

Caution! Never leave your child in the stroller while making adjustments! It can put your child in danger.

10.4. Adjustable backrest



Fig 10.3
Reclining back handle

The stroller's backrest can be adjusted to various positions varying from a seating position to a lying position. To change the position, move the backrest towards a seating position first. Then, pull the handle located at the back of the backrest and change the position slowly. After reaching the desired position, pull the handle down.

Caution! Each change of the backrest position requires an adjustment of the 5-point belts.

Caution! Never leave your child in the stroller while making adjustments! It can put your child in danger.

11. Basic equipment

11.1. Five-point belts



Fig 11.1.
Five-point belts

The Mewa stroller is standard equipped with professional five-point seat belts

- Five-point belts adjustment. Adjust the five-point belts to your child before attaching them. To adjust the length of the belts, move the adjusters on the braces up or down and move the protectors accordingly. To adjust the height of the braces to be attached to the backrest, move the braces from low to high, or vice versa.

- Use. In order to put on the child five-point belts, put the braces on the right and left of the child and then attach both braces to the belt buckle, which is fixed on the seat. In order to remove five-point belts from a child, press the button on the buckles and remove both braces from the buckle and slowly remove these braces from the right and left side of the child.

Attention! Each change in the angle of the backrest requires a readjustment of the five-point belt length.

11.2. Safety barrier

The Mewa stroller can be equipped with a safety barrier that prevents the sitting child from falling forward in certain indications..

- The safety barrier is mounted with the same fixings as the therapeutic table, so you cannot use both of these elements at the same time. To fastening the safety barrier to the frame slide it into the attachment tubes and listen for the 'click'. The clearly audible 'click' is an indication that the barrier has been securely attached to the seat. To disassemble the barrier pull both barrier latches on the frame and pull out the barrier at the same time.



Fig 11.2
Safety barrier

11.3. Adjustable canopy with window

The Mewa stroller can be equipped with an adjustable canopy with a window that helps protect the child from the sun or rain. The elements of the canopy are: canopy frame (front and rear arch), canopy opening adjustment mechanism (right and left), two snaps with red buttons fastening the canopy to the frame, canopy upholstery, and a mosquito net with a zipper.

- Assembly. The fastening of the adjustable canopy with the window are two plastic elements located on the outer sides of the upper part of the backrest. These attachments are used to mount the canopy with the window.

To install the adjustable canopy with the window, place the backrest upright, and then slide both canopy latches into the backrest mounts in parallel. Mount the canopy in such a way that the upholstery at the back/bottom of the canopy overlaps the back of the stroller. In

addition, make sure that the canopy adjustment element is secured with upholstery inside. This can be done by stretching the upholstery and placing the loop on the button in the corner of the canopy.

- Disassembly. To disassemble the adjustable canopy with window, press the two red buttons on the roof fixing latches and slide the latches out of the backrest mounts.

- Use. To change the size of the canopy, hold the front arch and the rear arch of the canopy frame and bring the arches closer or closer to each other to set the desired position. The canopy surface adjustment mechanism latches automatically. A special feature of the canopy is a window with a mosquito net for ventilation. The window opens and closes with a zipper.



Fig 11.3
Adjustable canopy

11.4. Shopping basket

The Mewa stroller can be equipped with a shopping basket which is used to transport small, ordinary items of everyday life, resistant to shocks, impacts, bruises, soaking, flooding and dirt.

It is mounted at the bottom of the stroller by wrapping it on the frame of the stroller and fastening it with a zipper.

The maximum load of the shopping basket is 3 kg.



Fig 11.4
Shopping basket

11.5. Abduction and stabilizing belts



Fig 11.5

Abduction and stabilizing belts

Abduction and stabilizing belts are used for additional pelvic stability and prevent the child from slipping out of the stroller. The correct way to use them is to fasten the child's thighs so that his pelvis adheres as close as possible to the backrest. This will ensure a healthier and more comfortable position of the child in the stroller.

- **Adjustment.** To adjust the length of the Abduction and stabilizing belts, you must move the adjusters on the front belt straps up or down.

- **Use.** In order to put the child on the Abduction and stabilizing belts, put the child on the installed belts, then wrap the child's leg and

fasten both belt buckles to both front straps. In order to remove the Abduction and stabilizing belts from the child, unfasten both belt buckles from the front straps.

11.6. Wedge

The Mewa stroller can be equipped with a wedge which in the case of specific indications provides discharge (inclination to the sides) of the child's lower limbs from his body. The wedge has a two-stage adjustment level. Fig.11.6



Fig 11.6
Wedge

11.7. Lining pads

The Mewa stroller is equipped with lining pads, securing the velcro on the trolley's backrest, if no additional equipment is installed. Fig 11.7.



Fig 11.7

Lining pads

11.8. Upholstery

The Mewa stroller is equipped with upholstery, which is completely removed for cleaning. To do this, detach all the elements fastening it to the stroller frame. How to clean the upholstery is described in point 17.8. Cleaning and care.

12. Additional equipment

12.1. Adjustable head supports - ↔ and ⇅

To adjust the head support, draw aside the seat upholstery on the front and back of the seat. Next, take off the head support knobs, which are situated on the back of the seat and adjust the head supports at the front of the seat. That accomplished, screw the knobs back in the appropriate place at the back of the seat. Finally, draw the upholstery at front and back of the seat. Fig 12.1



Fig 12.1
Side pads

12.2. Adjustable chest and loin pads - ↔ and ⇅

To adjust the chest and loin pads, draw aside the seat upholstery on the front (Fig 12.2.1) and back of the seat (Fig 12.2.2). Next, take off the chest and loin pads knobs, which are situated on the back of the seat (Fig 12.2.3) and adjust the chest and loin pads at the front of the seat. That accomplished, screw the knobs back in the appropriate place at the back of the seat. Finally, draw the upholstery at front and back of the seat. (Fig 12.2.4)



Fig 12.2.1

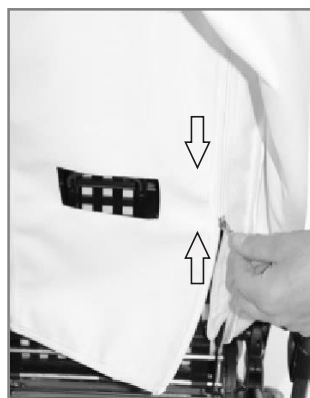


Fig 12.2.2



Fig 12.2.3



Fig 12.2.4

Unfasten the upholstery at both sides / Release the knobs / Fasten the upholstery



Fig 12.3
Therapeutic table

12.3. Therapeutic table

The Mewa stroller can be equipped with a non-slip therapeutic table with a maximum load of 3kg. – Fig.12.3.

The therapeutic table is mounted with the same fixings as the safety barrier so you cannot use both of these elements at the same time.

Slide the table into the attachment tubes and listen for the 'click'. To assemble the table, slide them onto the frame, latch. The table will not be properly installed until the latch engages the frame and snaps into place. This is evidenced by the clear "click" sound

coming from a well-connected therapeutic table mount. To disassemble the table, pull both latches on the frame and pull out the table at the same time.

12.4. **Seat minimizing kit**

The Mewa stroller is equipped with set of five inserts (Fig. 12.4). They are used to adjust the inside dimensions of the seat width. The upholstery of the seat minimizing kit is removable thanks to the zipper sewn in, which allows it to be replaced or cleaned.

- **Assembly.** The inserts are installed anywhere in the stroller, at any height adapting to the patient's indications. Installation is done with Velcro on the bottom, back and outer side of the insert.

- **Disassembly** is done by removing the element - pulling the velcro from the stroller upholstery.



Fig 12.4
Seat minimizing kit

12.5. **Travel bag**

The large and elegant travel bag is mounted on two handles on the inside of the handle to push the stroller (see Fig. 12.5). The maximum bag load is 2kg.

12.6. **Sun umbrella**

The Mewa stroller can be equipped with a Sun umbrella (see Fig. 12.6) that helps protect the child from the sun.

12.7. **Winter gloves**

Two winter gloves that increase the comfort of a child's assistant during colder weather can be put on the handle to push the stroller (see Fig.12.7). First you need to put on two material clamps and fasten them with velcro on the handle to push the stroller, then gloves are put on these material clamps and fastened with snaps.



Fig 12.5
Travel bag



Fig 12.6
Sun umbrella



Fig 12.7
Winter gloves

12.8. Winter sleeping bag

The Mewa stroller can be equipped with a Winter sleeping bag, increasing the child's comfort during colder weather (Fig. 12.8). This optional equipment requires installation to other optional equipment - a barrier. To mount the Winter sleeping bag, place it on the footrest from below, and transfer it through the safety barrier from above and fasten it with snaps to the frame of the stroller.

12.9. Mosquito net

The Mewa stroller can be equipped with a mosquito net (see Fig. 12.9), which prevents insects from coming into direct contact with the child. This additional accessory requires fitting to other optional equipment. In order to mount the mosquito net, put it from above on the entire adjustable canopy with window, and then cover the mosquito net seat and footrest.

12.10. Raincover

The Mewa stroller can be equipped with a rain cover that helps protect the child from moisture. This optional equipment requires mounting to other optional equipment. To mount the rain cover, put it on top of the entire adjustable canopy with window and then cover the seat and footrest with the cover (see Fig. 12.10).



Fig 12.8
Winter sleeping bag



Fig 12.9
Mosquito net



Fig 12.10
Raincover

13. Technical specification of the Mewa stroller

Manufacturer	MyWam Kupiec, Bartold, Angres spółka jawna
Address	ul. Szczecińska 10, 41-516 Chorzów, Poland
Type	Special children's stroller
Model	Mewa

PARAMETERS	pomiar
stroller length	111 cm
stroller width	65 cm
stroller height	120 cm
total basic mass	23 kg
effective seat depth	26,5*/29
seat width	28*/33
back height	62 cm
backrest width	31 cm
distance between footrest and seat	19-29 (every 2,5 cm)
folded stroller length for transport	frame 84cm / seat 88cm
folded stroller width for transport	frame 64 cm/ seat 45cm
height of the folded stroller for transport	frame 26 cm/ seat 43cm
basket load capacity	3 kg
bag capacity	2 kg
maximum load of the stroller	35 kg

The manufacturer reserves the right to make changes to technical parameters. Measurement values tolerance +/- 20mm, +/- 1.5kg, +/- 1°.

The stroller Mewa has been positively tested in terms of standards:

PN-EN 12183:2014	PN-EN 12182:2012	ISO 7176-1:2014
ISO 7176-3:2012	ISO 7176-5:2008	ISO 7176-7:2001
ISO 7176-8:2014	PN-EN 1021-1:2007	PN-EN 980:2010
PN-EN 1041:2010	PN-EN ISO 14971:2011	

13.1. User parameters

max. patient weight	35 kg
patient height	90-125 cm

13.2. Material composition of the stroller

Powder coated steel and aluminium, fabrics 100% polyester and plastic (PP, PCV, ABS, PA, PE, PU).

14. Delivery and tools

The Mewa stroller is delivered in a cardboard box and pre-assembled. The original packaging should contain the following items:

- Folded stroller frame with seat, backrest and wheels
- Optional equipment as ordered

- User manual with warranty card

Before use, check and make sure that the carton contains all the elements and none of them has been damaged, e.g. during transport. In the event that it is found that in the carton received is missing or damaged any component, contact the Seller and the Manufacturer in order to receive the missing elements or replace the product with a new, fully complete and functional.

The multifunctionality and extensive equipment of the Mewa stroller allow you to tailor the stroller to your personal requirements. To do this, read the entire manual carefully before use.

No tools are required to assemble and disassemble the Mewa stroller.

A Phillips screwdriver is required to adjust the height of the footrest. An Allen key is needed to adjust the backrest height and seat depth.

15. Assembly and disassembly of the stroller

- **Stroller assembly.** The product pre-packed in the original carton is already pre-assembled. Complete the installation by completing the following steps:
 - a) mount the rear wheels on the axle
 - b) unfold the stroller frame with a footrest adjusting the parameters:
 - enable/disable lock of the direction of travel of the front wheels
 - adjust the height and angle of the stroller push handle
 - adjust the seat tilt angle
 - set the tilt angle of the backrest
 - adjust the length of the five-point belts
 - adjust the height of the footrest and the tilt angle of the footrest
 - c) install appropriate accessories and adjust parameters, for example:
 - adjust the usable width and depth of the seat (set minimizing kit)
 - install a safety barrier or a therapeutic table
 - adjust the deflection and stabilization belts
 - adjust the canopy with the window
 - install and adjust the wedge
 - put on a travel bag
 - put on a shopping basket
 - put on selected items of equipment dedicated to the seasons: winter sleeping bag, winter gloves, mosquito net, raincover, umbrella
- **Dismantling the stroller.** In order to dismantle the stroller, it is necessary to:
 - a) remove the optional equipment
 - b) unclip the seat from the frame of the stroller
 - c) lock the front wheels, fold the footrest and stroller frame
 - d) unclip the rear and front wheels
 - e) fold in the stroller closing it with the transport lock

16. Preparation for use

Before using the stroller:

- it is necessary to ensure that the child's footwear is comfortable and stabilizes the ankle well.
- make sure that the child is not placed in the stroller without clothing covering the body in places of contact with the upholstery of the stroller. This is to prevent chafing of the child's body and hygiene of the child.
- For children who have problems with urine retention, put a special oilcloth on the seat of the stroller, which can be purchased in pharmacies or medical stores.
- When using any additional aids such as venflons, drains, bags, catheters, etc., make sure that they are well secured against falling out, pulling out, clamping and do not interfere when using the stroller.

In order to start using the stroller properly, all the operations listed below must be carried out:

- Check the good technical condition of the stroller and all its components, and if the technical condition of the stroller or its individual component is bad or the user has any doubts about it, then the stroller must not be used. Please contact the Seller and the Manufacturer.
- Place the stroller on a flat and secure surface.
- Unfold the stroller correctly
- Engage the central brake.
- Make all necessary adjustments to the stroller
- Prepare for five-point belts, straps on the footrest and Abduction and stabilizing belts, if these accessories are available.
- Place the child in the stroller and make sure that all adjustments previously made are appropriate and safe for the child.
- Attach the child to five-point belts and straps on the footrest.
- Put on the child other necessary stroller equipment, if available.

17. Using the stroller

17.1. Carrying the stroller

The best way to carry a stroller is to use it with its wheels. In cases where this is not possible (e.g. to transport the stroller in a car), for safety reasons, the stroller can be lifted and moved after dismantling the additional equipment and the front and rear wheels, in order to reduce the load. When carrying or lifting the stroller frame with the wheels mounted, fold the stroller frame and then hold the stroller frame only by the frame's fixed elements. For safety reasons, do not lift the entire stroller, i.e. the frame with the seat and backrest, wheels, footrest and mounted accessories, as this may damage your health. Do not lift the stroller by its footrest or safety barrier, therapy table or any other accessories, as this may damage the product and endanger the safety of others.

17.2. Getting on and off the stroller

When getting on and off the stroller, place the stroller as close as possible to a stable and safe place where the user will change. After making sure that the central brake is on and the footrest is raised to avoid standing on it, the user can change.

17.3. Correct position in a stroller

For safety and health reasons, it is recommended that the user have the pelvis and back as close as possible to the backrest as much as possible to support the spine evenly (optimal support for all parts of the spine so that there is no free space).

The thighs should be supported on the largest possible surface, along their entire length, which is used to adjust the depth of the stroller seat, and the feet should be supported on the footrest plate which prevents slipping, maintains optimal leg muscle tension and carries some of the weight from the upper body.

17.4. Riding in a stroller

When using the stroller, always hold the handle to push the stroller to prevent the wheel from rolling slowly. Release the central brake while holding the push handle of the stroller and push the stroller in a safe direction.

17.5. Moving on slopes

Always position the stroller forwards when climbing or driveways. Limit the speed on sloping surfaces and move as slowly as possible. The stroller assistant (parent, guardian) must have sufficient strength to control the stability of the stroller.

17.6. Overcoming thresholds or curbs

When overcoming small obstacles (small thresholds, small curbs, etc.), block the tilting foot with the foot and the handle to push the stroller downwards and towards you to raise the front wheels of the stroller and carefully drive the front wheels over the obstacle. Then slowly pull the push handle upwards away from the ground to raise the rear wheels of the stroller and carefully drive the rear wheels onto the obstacle.

17.7. Transport in a car and other vehicles

To transport the Mewa stroller in the car, remove all accessories from the stroller and store them in a safe place. Then dismantle the stroller (see chapter 15) and place it in the trunk.

If the dismantled stroller is not in the trunk of the car separated from the passenger compartment, then the stroller frame must be securely attached to the vehicle. For this purpose, unused seat belts may be used. In order to transport a stroller in other vehicles, such as trains, trams, buses, planes, ships, etc., please ask the persons responsible for transport in these means of transport about the conditions of safe travel to the destination with a stroller and the appropriate requirements.

In order to transport the stroller on other vehicles, remove all accessories from the stroller and store them in a safe place. Then dismantle the stroller (see chapter 15) and place it in the baggage storage area of the given vehicle. If the dismantled stroller is not in the trunk of the vehicle separated from the passenger compartment, the stroller frame must be securely attached to the vehicle.

The Mewa stroller should not be used as a seat to transport a child in a vehicle.

17.8. Cleaning and care

Small leaks and stains can be removed from the surface with a kitchen sponge and using clean, warm water.

Stains should be removed immediately after they appear on the surface of the stroller. Do not rub - rubbing will spread the stain to a larger surface and may distort the texture of the fabric.

The treatment of water stains such as stains from food, drinks and vegetable oils is as follows:

- Dry the stain by pressing a white paper towel, handkerchief or absorbent clean cloth over its entire surface.
- Soak a small piece of cloth with a weak solution of detergent and warm water (or cold water in the event of blood staining) and rinse the stain by removing it from the outside to the inside to prevent it from growing. Do not use a too wet cloth, as excessive amount of liquid may cause ripples after drying or enlargement of the stain and damage of the upholstery or foam inside.
- It may be necessary to repeat the above operation depending on the type of stain and the amount of stained material. The above operation should be continued until the stain is completely removed, but always remember to remove excess moisture from the stained material before rinsing again.
- Then rinse the material using only lukewarm water without detergent, using a clean sponge or a clean white cloth. As above, rinse the material in place after the stain in the right direction, i.e. from the outside in.

Disinfectants that do not contain organic solvents can be used in accordance with the manufacturer's instructions.

Do not dip the fabric. Make sure the product is always dry before use.

Do not use organic solvents or dry cleaning agents.

For disinfection, you can use a small amount of ordinary household bleach, diluted with water with a maximum concentration of 3% bleach.

ATTENTION! The manufacturer is not liable for damages resulting from the use of inappropriate cleaning products.

18. Storage of the stroller

The stroller should be stored in a clean, dry and well ventilated room, preferably at 0-35°C.

18.1. Preparation of the stroller for long-term storage

In order to prepare the Yeti stroller for long-term storage, secure the product and perform all service operations in the "Periodic Inspections" chapter as if they had occurred.

The stroller should be stored in a clean, dry and well ventilated room.

18.2. Corrosion protection

To avoid corrosive damage to the stroller, it should not be used or left in wet and damp rooms.

When the truck is wet, wipe all wet parts with a soft cloth and dry with a light oil preservative. This does not apply to the upholstery of the stroller, the provision of cleanliness and maintenance of which is described in detail in the Cleaning and care chapter.

As a standard, the stroller is covered with the highest quality coatings. In the event of damage to the varnish, protect the damage by applying a touch-up compound using the varnish proposed by the varnish Manufacturer.

Regular cleaning and lubrication of the stroller mechanisms with light oil guarantees long and trouble-free use of the stroller.

19. Re-use

The product is suitable for reuse.

All reusable products are subject to a particular burden.

Therefore, their features and performance may deviate from the new product to such an extent that the safety of patients and third parties during their use is at risk.

The product with proper care and maintenance retains its quality for a period of 3 years and this is the expected period of its use.

For secondary use, please:

- Remove all upholstered parts and replace them with new ones.
- The stroller frame with wheels and steering elements should be thoroughly cleaned and disinfected with commercially available cleaning and disinfecting agents.
- The product must then be checked by an Authorized Service Center (see Chapter 25) or the Manufacturer for wear and damage.
- All worn parts and components that do not match or are unsuitable for the user should be replaced with new ones.
- The product prepared in this way can be put back into use.

20. Service and periodic inspections

For safe and trouble-free use of the Yeti stroller, it must be regularly and thoroughly inspected and serviced as indicated in the table below.

In a situation where the stroller user is not able to carry out the activities himself, he should return the stroller to a specialized stroller service point or directly to the Manufacturer, who will properly review and will perform necessary service activities. These activities may be charged extra.

All repairs should be carried out in specialized, authorized Services (see point 25) or directly at the Manufacturer.

Only original spare parts that are supplied or recommended by the Manufacturer may be used in the stroller.

Service should be performed:

- **Each time before using the stroller**
 - Check the efficiency of the central brake - After activating the central brake, both rear wheels of the stroller should be locked and cannot rotate.
 - Check the attachment of structural elements - It is necessary to check all mechanisms of the stroller and its movable elements, in particular, pay attention to: the correct mounting of the therapeutic seat with the backrest to the stroller frame, the correct mounting of all four wheels, the correct functioning of the five-point seat belts, properly locked seat angle adjustment mechanism, properly locked backrest angle adjustment mechanism and footrest correct mounting.
 - Check the condition of the tires - Check that the tires on all four wheels are not damaged or worn.
 - Inspection of solid frame parts, seat with backrest and footrest frame - Visually inspect the condition of the structural profiles. These elements must not be damaged or cracked.
- **Every one month**

- Cleaning the stroller - Clean any moving parts, bearings, adjustment mechanisms with a soft cloth and then lubricate them with a few drops of light oil. Be careful not to soil the upholstery of the stroller or your own clothes.
- Check all bolted connections for looseness - Carefully check that all bolts and knobs in the stroller are properly tightened.
- **Every six months**
 - Check the cleanliness of the upholstery and inspect the general condition of the stroller - Do not allow the stroller to be neglected as described in the chapter „Cleaning and care“.
 - Checking the general condition of the stroller. The stroller must not be neglected





21. Disposal

Disposal and recycling: the used stroller should be returned to a recycling center in accordance with the regulations in force in the given country.

22. Compliance with requirements for medical devices

The Yeti stroller meets the essential requirements of Council of Europe Directive 93/42 EEC amended by 2007/47/EC for medical devices.

The Yeti stroller, in accordance with Directive 93/42 EEC, is a non-invasive, inactive medical device of class I according to rule 1.

	<p>MYWAM Kupiec, Bartold, Angres spółka jawna ul. Szczecińska 10, 41-516 Chorzów, Polska tel. 0048 32 733 11 31 mywam@mywam.pl</p>													
<p>DEKLARACJA ZGODNOŚCI WE EC DECLARATION OF CONFORMITY</p>														
<p><u>Wytwórca:</u> <u>Manufacturer:</u></p>	<p>MYWAM Kupiec, Bartold, Angres spółka jawna ul. Szczecińska 10, 41-516 Chorzów, Poland</p>													
<p><u>Produkt:</u> <u>Product:</u></p>	<p>Wózek inwalidzki specjalny MEWA Specialist Stroller MEWA</p>													
<p><i>Deklaruję z pełną odpowiedzialnością, że wózek inwalidzki specjalny Mewa, będący wyrobem medycznym klasy I, jest zgodny z wymaganiami dyrektywy Rady Europy o wyrobach medycznych - 93/42/EEC z 14 czerwca 1993 roku zmienionej przez 2007/47/EEC z 21 marca 2010 roku, przyjętej w Polsce na mocy ustawy o wyrobach medycznych z dnia 20 maja 2010 r. (Dz. U. nr 107, poz. 679), oraz zgodny z normami zharmonizowanymi:</i></p>														
<p><i>I solely declare that the specialist Stroller Mewa is a class I device, and it is in conformity with the requirements of EU Council Directive for Medical Devices 93/42/EEC of 14 June 1993 amended by 2007/47/EEC of 21 March 2010, implemented in Poland under the Medical Devices Act of 20 May 2010 (Dz. U. nr 107, poz. 679), and is in conformity with harmonized standards:</i></p>														
<table border="1"> <tr> <td>PN-EN 12183:2014</td> <td>PN-EN 12182:2012</td> <td>ISO 7176-1:2014</td> </tr> <tr> <td>ISO 7176-3:2012</td> <td>ISO 7176-5:2008</td> <td>ISO 7176-7:2001</td> </tr> <tr> <td>ISO 7176-8:2014</td> <td>PN-EN 1021-1:2007</td> <td>PN-EN 980:2010</td> </tr> <tr> <td>PN-EN 1041:2010</td> <td>PN-EN ISO 14971:2011</td> <td></td> </tr> </table>	PN-EN 12183:2014	PN-EN 12182:2012	ISO 7176-1:2014	ISO 7176-3:2012	ISO 7176-5:2008	ISO 7176-7:2001	ISO 7176-8:2014	PN-EN 1021-1:2007	PN-EN 980:2010	PN-EN 1041:2010	PN-EN ISO 14971:2011			
PN-EN 12183:2014	PN-EN 12182:2012	ISO 7176-1:2014												
ISO 7176-3:2012	ISO 7176-5:2008	ISO 7176-7:2001												
ISO 7176-8:2014	PN-EN 1021-1:2007	PN-EN 980:2010												
PN-EN 1041:2010	PN-EN ISO 14971:2011													
<table border="1"> <tr> <td>ISO 7176-3:2012</td> <td>ISO 7176-5:2008</td> <td>ISO 7176-7:2001</td> </tr> <tr> <td>ISO 7176-8:2014</td> <td>PN-EN 1021-1:2007</td> <td>PN-EN 980:2010</td> </tr> <tr> <td>PN-EN 1041:2010</td> <td>PN-EN ISO 14971:2011</td> <td></td> </tr> </table>	ISO 7176-3:2012	ISO 7176-5:2008	ISO 7176-7:2001	ISO 7176-8:2014	PN-EN 1021-1:2007	PN-EN 980:2010	PN-EN 1041:2010	PN-EN ISO 14971:2011						
ISO 7176-3:2012	ISO 7176-5:2008	ISO 7176-7:2001												
ISO 7176-8:2014	PN-EN 1021-1:2007	PN-EN 980:2010												
PN-EN 1041:2010	PN-EN ISO 14971:2011													
<table border="1"> <tr> <td>ISO 7176-8:2014</td> <td>PN-EN 1021-1:2007</td> <td>PN-EN 980:2010</td> </tr> <tr> <td>PN-EN 1041:2010</td> <td>PN-EN ISO 14971:2011</td> <td></td> </tr> </table>	ISO 7176-8:2014	PN-EN 1021-1:2007	PN-EN 980:2010	PN-EN 1041:2010	PN-EN ISO 14971:2011									
ISO 7176-8:2014	PN-EN 1021-1:2007	PN-EN 980:2010												
PN-EN 1041:2010	PN-EN ISO 14971:2011													
<table border="1"> <tr> <td>PN-EN 1041:2010</td> <td>PN-EN ISO 14971:2011</td> <td></td> </tr> </table>	PN-EN 1041:2010	PN-EN ISO 14971:2011												
PN-EN 1041:2010	PN-EN ISO 14971:2011													
<div style="display: flex; justify-content: space-between; align-items: center;"> <div data-bbox="389 1704 509 1794">  </div> <div data-bbox="1050 1756 1246 1827"> <p>Małgorzata Uchman <i>[Signature]</i> KIEROWNIK DZIAŁU SPRZEDAŻY</p> </div> </div>														
<p>Chorzów, 05 listopad 2019r./05 November 2019</p>														
<p>Data/Date</p>	<p>Podpis/Signature</p>													
<div style="display: flex; justify-content: center; align-items: center;">  <div data-bbox="1050 1912 1278 2018"> <p>MYWAM Kupiec, Bartold, Angres Spółka Jawna 41-516 Chorzów, ul. Szczecińska 10 www.MyWam.pl tel. 32 733 11 31, email: mywam@mywam.pl NIP 627 273 67 93, REGON 243428050 KRS 0000487471</p> </div> </div>														

23. Notes

24. List of authorized service points

POLAND

MyWam Kupiec, Bartold, Angres sp. j.

ul. Szczecińska 10, 41-516 Chorzów, PL

tel: +48 32 733 11 31

www.mywam.eu

mywam@mywam.pl

25. Warranty Card

Series and number:

Date of purchase: - -

Seller stamp and signature

26. Warranty conditions

The manufacturer grants the buyer a warranty for the purchased stroller, for a period of 24 months from the date of purchase.

Defects revealed during this period will be removed free of charge.

The Buyer submitting a complaint on the stroller should deliver the stroller to the Seller or to the Manufacturer in the original factory packaging, together with the elements of the stroller protected against damage.

Damage caused by insufficient securing of the stroller shipment and its components is not subject to warranty repairs.

The complainant may not demand gratification if he improperly marks or incorrectly insures the shipment of the stroller and its components.

The cost of shipping the complained stroller to the Seller or to the Manufacturer shall be borne by the buyer complaining the stroller.

The warranty period is extended by the time needed to repair the stroller.

The activities provided for in the manuals, which the user is obliged to perform on their own and at their own expense, cannot be performed as part of a warranty repair.

The buyer has the right to replace the complained stroller with a new stroller if the same damages or defects that have already been removed twice and which prevent the stroller from being used are still present.

The stroller should be used in accordance with these operating manual.

Unauthorized use of the stroller, inconsistent with these operating manual or inconsistent with the principles of correct operation will void the warranty.

Making any structural changes to the stroller will void the warranty.

Damage to the stroller caused by: improper operation, improper maintenance, incorrect storage of the stroller, mechanical damage as well as natural abrasion and wear of parts (e.g. upholstery, wheels) are not covered by the warranty. They are repaired at the user's expense.

The manufacturer warns against performing warranty repairs at unauthorized service facilities. Finding such repair or unauthorized alteration of the stroller will void the warranty.

The warranty card is valid provided that it has: the date of sale entered, confirmed by the Seller's company seal, the Seller's signature, type and serial number of the stroller.

The warranty card is valid as long as it does not contain deletions and corrections.

The stroller manufacturer runs its paid post-warranty service.

In matters not covered by this warranty card, the provisions of the Civil Code shall apply.

All matters related to the warranty and complaint of the stroller should be directed to the Seller or directly to the Manufacturer.