



CEMD MEDICAL DEVICE - CLASS I

SPECIAL NEEDS STROLLER FOR CHILDREN

GENERAL INFORMATION

The YETI special needs stroller is designed for larger/older the children who are unable to walk or sit in a correct position as a result of injuries, diseases or motor dysfunctions. The stroller is suitable for users at a height of 130 - 175 cm and a maximum weight of 80 kg.

Its heavy-duty constructions is equipped with the "I Grow With You" system (smooth extension of the backrest height and seat depth as you child grows) and adjustable tension system which allows to adjust the proper upholstery tension both on the seat and the backrest guarantees you a support for many years. Moreover, **YETIs** solid suspension along with a fast folding mechanism, tilted seat, reclining backrest and a wide range of accessories significantly simplify the everyday use irrespective of the weather conditions and time of the year. Modern and practical design is his additional asset.

INDICATIONS

paralysis and paresis of lower limbs of various origins, paralysis and paresis of the trunk, weakening of postural muscles responsible for maintaining the correct body position, cerebral palsy, myelomeningocele, muscular dystrophies.


CONTRAINDICATIONS

In case of any contraindications to sitting position &/or passive verticalization stroller adjustments should be consulted with physiotherapist or physician-in-charge.

NEW

FEATURES / SAFETY AND COMFORT

TILT IN SPACE SEAT simple mechanism RECLINING BACKREST wide range of adjustments ADJUSTABLE AND ANGLE TILT FOOTREST simple and quick adjustment

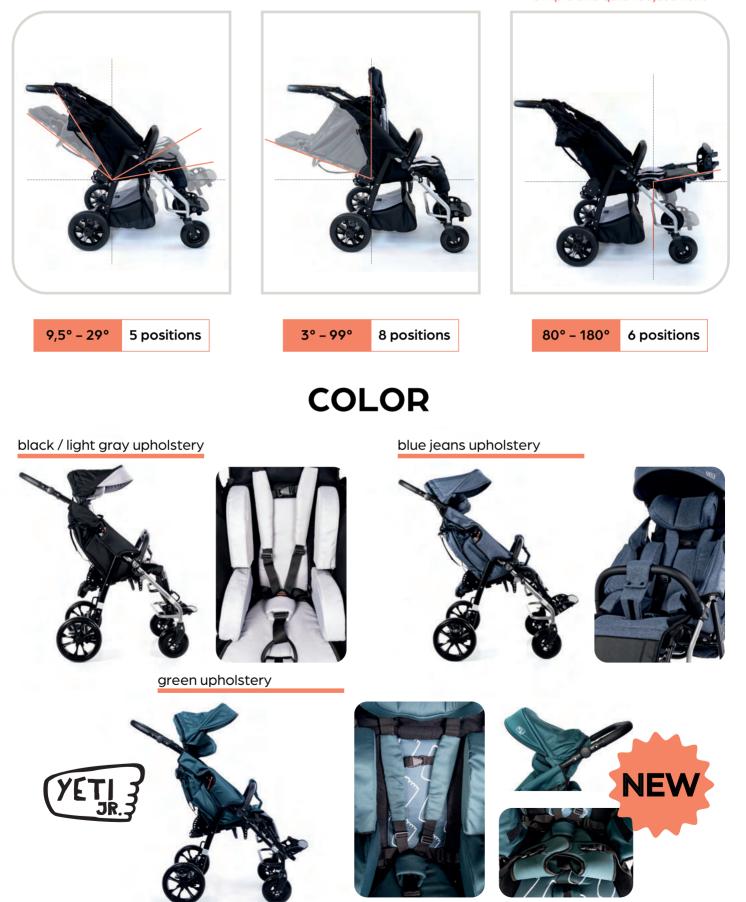

WINDSTOPPER

multifunctional upholstery

FOLDING LIGHTWEIGHT FRAME

ease of storage and transportation

DRIVING SUPPORT comfort for the bumpiest journeys



Yeti is equipped with a multifunctional windstopper that not only comes in handy during windy conditions but also protects against unnecessary arm movement and limits visual stimulation.



The one-piece construction **Yeti** stroller allows for a quick and easy folding which facilitates daily use, storage and transport.



The adjustable tension system uses 14 independent velcro straps (7 on the seat and 7 on the backrest), to customize **support and stiffness** to the child's needs by tightening or loosening specific points.

FUNCTIONAL ACCESSORIES IN PRACTICE





TECHNICAL PARAMETERS OF THE STROLLER

seat width: 28*/43 cm effective seat depth: 31 – 45 cm backrest height: 65 – 85 cm backrest width: 25*/42 cm distance between footplate and seat: 28,5 – 43,5 cm, 7 positions *values are given with inserts total weight: **23 – 26 kg** stroller length: **114,5 – 146 cm** stroller height: **117 – 133 cm** stroller width: **67 cm**

IMPORTANT NOTE! The choice of a medical device appropriate for a specific disease entity and condition of the patient must be made with the help of a specialist doctor or physiotherapist. Before use, read the operating manual and the label or consult your doctor, as this medical device may not be suitable for you.