

# mouse

## SPECIAL NEEDS STROLLER FOR CHILDREN



  MEDICAL DEVICE - CLASS I

## GENERAL INFORMATION

The **MOUSE special needs stroller** is designed for the youngest users who are unable to walk or sit in a correct position as a result of injuries, diseases or motor dysfunctions. The stroller is suitable for children at a height of **80 – 110 cm** and a maximum weight of **30 kg**.

Thanks to the backrest tilt and recline feature as well as adjustable torso and head pads you may easily adjust **MOUSE** to the required position and ensure the user with optimal protection and stability.

The Mouse stroller features rear-wheel suspension and off-road wheels, making it much easier to travel over uneven terrain.

## INDICATIONS

Paralysis and paresis of lower limbs of various origins, paralysis and paresis of the trunk, weakening of postural muscles responsible for maintaining the correct body position, cerebral palsy, myelomeningocele, muscular dystrophies.

## CONTRINDICATION

In case of any contraindications to sitting position &/or passive verticalization stroller adjustments should be consulted with physiotherapist or physician-in-charge.

**THIS IS MEDICAL DEVICE.**

**USE IT IN ACCORDANCE WITH THE INSTRUCTIONS FOR USE OR THE LABEL.**

# FEATURES / SAFETY AND COMFORT

## TILT IN SPACE SEAT



8° – 23°

4 positions

## RECLINING BACKREST



16° – 88°

4 positions

## ADJUSTABLE AND ANGLE TILT FOOTREST



84° – 173°

6 positions

## SEAT MINIMIZING KIT reducing insert



The **seat minimizing kit** in **Mouse** is a comfortable, one-piece insert that allows you to reduce both width and depth of the seat.

## STABILIZATION SYSTEM adjustable pads



The **head and torso support** consists of two pairs of independent pads that are height, width and angle adjustable. Its wide range of adjustments and simplicity gives you a fast and effective adaptation to the user's individual needs.

## LYING POSITION



With a fully adjustable seat and footrest you may quickly alter the stroller into cozy bed, perfect for a short rest or a nap.



**REVERSIBLE SEAT**

arrange your travel easily

**FOLDING LIGHTWEIGHT FRAME**ease of storage, transportation  
and dismantling**DRIVING SUPPORT**

comfort for the bumpiest journeys



The design of the stroller allows you to easily rotate the seat opposite to the driving direction. Thanks to this feature you can keep an eye on your child all the time.



The construction of main elements allows you to dismantle **Mouse** into 3 separate units (frame, seat and wheels).



The stroller is equipped with **a shock-absorbing suspension** on all four wheels and **an adjustable rear suspension** in the rear wheels with option to lock, **central brake** and **a tilt-foot**.

**FUNCTIONAL ACCESSORIES IN PRACTICE**


1 Canopy with window



2 Set minimizing kit



3 Torso/head pads

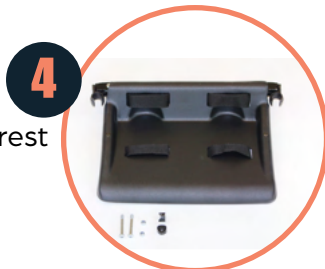

4 Adjustable footrest



5 Swivel front castor wheels with option to lock



6 Quick release mount on all 4 wheels



# EQUIPMENT



SET MINIMIZING KIT

HEAD/TORSO  
ADJUSTABLE PADS

ABDUCTION BELTS



SAFETY BARRIER



FIVE-POINTS BELTS

ADJUSTABLE  
ABDUCTION BLOCK

THERAPY TRAY

CANOPY WITH  
WINDOW

WINTER GLOVES



TRAVEL BAG



SHOPPING BASKET



SUN UMBRELLA



RAIN COVER



SLEEPING BAG



MOSQUITO NET

## TECHNICAL PARAMETERS OF THE STROLLER

seat width: 24\*/32 cm

effective seat depth: 22\*/26 cm

backrest height: 52 cm

backrest width: 30 cm

distance between footplate and seat: 14,5 – 23 cm, 4 positions

\*values are given with inserts

total weight: 17,2 – 23,2 kg

stroller length: 95 – 124 cm

stroller height: 122 – 134 cm

stroller width: 60 cm

**IMPORTANT NOTE!** The choice of a medical device appropriate for a specific disease entity and condition of the patient must be made with the help of a specialist doctor or physiotherapist. Before use, read the operating manual and the label or consult your doctor, as this medical device may not be suitable for you.